


*cryo***SPHERE[®]+**




Precision.
Protection.
Performance.

AtriCure

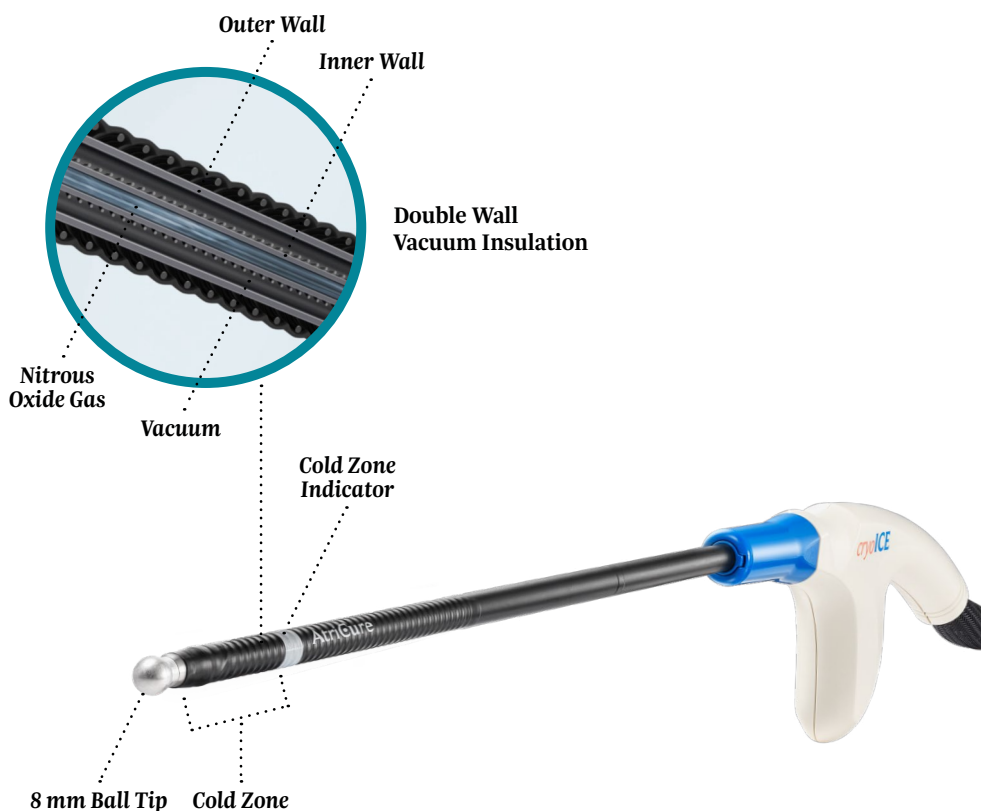

Cryoablation Probes

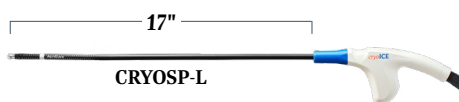

cryoICE BOX® - V6


ACM (US)

cryoICE Cart & Upgrade Kits

ASCS (cart with double tank holder)

20 lbs N₂O tanks not provided by AtriCure

Cryo Nerve Block: A non-opioid solution for post-operative pain

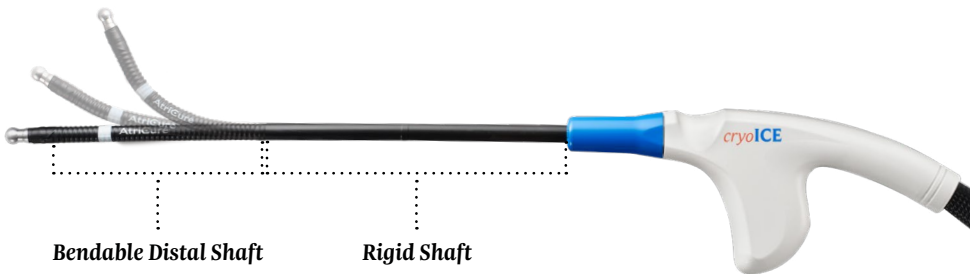

Features

Double Wall Vacuum Insulation	Designed to reduce risk of injury to non-target tissues that may contact the insulated shaft
Time To Therapeutic Temperature	Optimized time* to reach therapeutic temperature
Cold Zone Indicator	Provides visual reference that is approximately 4 cm* from the distal tip of the probe shaft**
8 mm Ball Tip	8 mm ball tip specifically engineered for cryoablation of peripheral nerves Purposefully shaped to reach target tissue from various approaches/angles
Active Defrost	Allows the probe to be easily removed from the cryoablation site
Rigid Shaft	Stainless steel shaft provides rigidity and tactile feedback for the user
Bendable Distal Shaft	Distal shaft bendable up to 140 degrees

*Reference on file at AtriCure

**Shaft distal of the Cold Zone Indicator can cause unintended cryoadhesion if touching unintended structures

U.S. Indications:

For Adult Patients AtriCure's cryoICE cryoSPHERE® + cryoablation probes are intended for use to temporarily block pain by ablating peripheral nerves performed by freezing target tissues, creating an inflammatory response (cryonecrosis). **For Adolescent Patients** The cryoICE cryoSPHERE+ cryoablation probes are intended for use to temporarily block pain by ablating intercostal nerves under direct visualization¹ in adolescent patients of at least 12 years of age.

¹Direct visualization, in this context, requires that the surgeon is able to see the targeted tissue for cryoablation directly or with assistance from a camera, endoscope or other similar optical technology.

ATRICURE, INC.

7555 Innovation Way
Mason, OH 45040 USA

+1 (513) 755-4100

+1 (866) 349-2342

CustomerService@AtriCure.com

www.AtriCure.com

PM-US-3132C-1027-G

AtriCure